

Dedicated in Ophthalmology

Precise Measurement with both Immersion & Contact mode
Touch Screen, Extremely Easy Use
Portable & Sleek Design
Auto / Manual Modes
Auto Gain Control
Built-in Thermal Printout



SPECIFICATIONS:

A-Scan (for MD-1000A & MD-1000 A/P)

- Probe: 10MHz with Fixation Red Light
- Total Gain: 100dB with an adjustable range of 0~50dB
- Biometry Accuracy: $\leq \pm 0.05\text{mm}$
- Resolution: 0.01mm
- Measuring Range: 15~40mm
- Measuring Mode: Contact, Immersion mode
- Measuring Parameters: Anterior Chamber Depth, Lens Thickness, Vitreous Length and Axial Length.
- 5 Different Measuring Modes: Automatic (for Normal, Cataract, Aphakic and Special), and Manual
- 8 Groups of Readings for Average, SD
- IOL Formula: SRK-T, SRK-II, BINKHOST- II, HOLLADAY, HOFFER-Q and HAIGIS

Standard Configuration:


- Main Unit:**
 - MD-1000A A-Biometer only
 - MD-1000P Pachymeter only
 - MD-1000A/P Biometer
- Probe:**
 - 10MHz A probe (MD-1000A, MD-1000A/P)
 - 20MHz P Probe (MD-1000P, MD-1000A/P)
- Footswitch
- Test Object
- PC Suite
- AC Adapter

Pachymeter (for MD-1000P & MD-1000 A/P)

- Probe Frequency: 15~20MHz
- Display Resolution: $1\mu\text{m}$
- Biometry Accuracy: $\leq \pm 5\mu\text{m}$
- Measuring Scope: 230~1200 μm
- Multiple Corneal Maps with Graphical Display

Others

- Power Supply: AC 100 ~ 240V, 50/60Hz, 50VA
- Dimension: 337mm x 177mm x 155mm (L x W x H)
- Weight: 1.7Kg

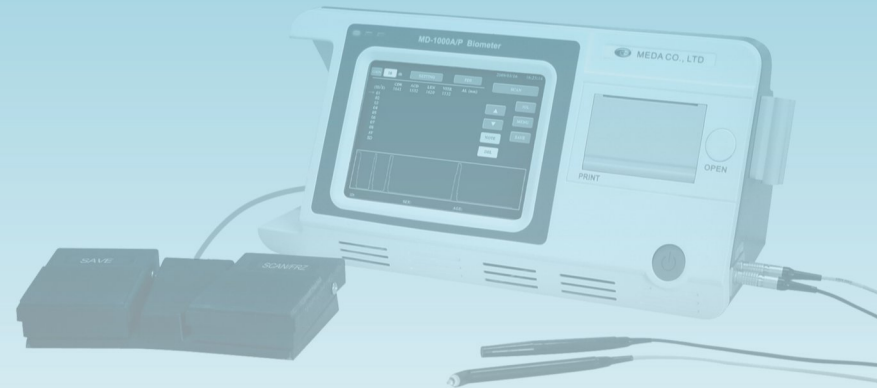

Room D,F3,Building C2,
Xinmao Science Skill Park,
Huayuan Industry
Development Area,
Tianjin 300384, China

Phone: +86-22-83713808 ext.0
+86-22-83713881
Fax: +86-22-83713880
www.MEDA.com.cn
E-mail: export@meda.com.cn

DISTRIBUTOR:



MD-1000A
MD-1000P
MD-1000A/P
Biometer for Ophthalmology






Precise & Accurate!

OPTIONAL



A Biometer Only

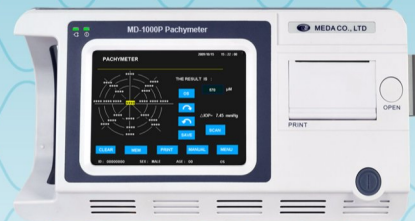


Pachymeter



A-Biometer / Pachymeter

A-SCAN

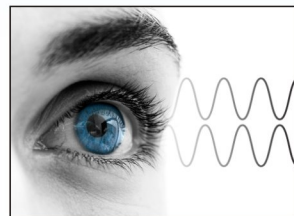

Precise and accurate

A scan measurement is ensured by MEDA's dedicated and matured technology in ophthalmic field.




Automatic Reading

System automatically measures 8 readings for one group. Users may manually adjust reading results during the automatic measuring.




Splendid Measuring Mode

Contact & Immersion measuring mode. Automatic reading with 4 different modes: Normal, Cataract, Aphakic and Special. Manual measuring available.

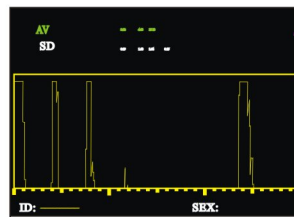


S.D. Function

MD-1000A (A/P) provides S.D.(Standard Deviation) to help evaluate reading reliability.

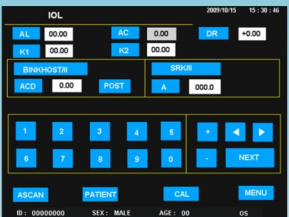

IOL Calculation


6 IOL formulae for IOL calculation. Instant formula switch. Automatically import Axial Length. Parameters input with Touch Screen.

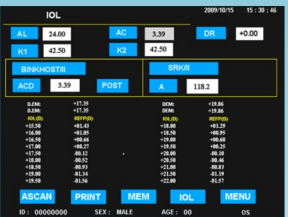

IOL Calculation


Dual formulae display for result comparison. Easy access to database. One button for instant printout.

PACHYETER

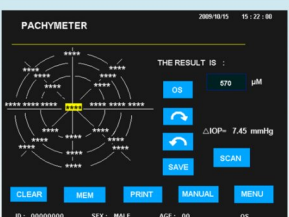


Accurate Measurement

Automatic reading for cornea thickness. Multiple or single point reading. Multiple measurements for single point cornea thickness. Average reading sharply increases accuracy.

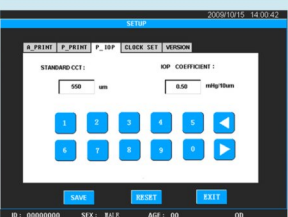


IOP (Intraocular Pressure) Adjustment

Provides reference for tonometer measurement. Parameters adjustable according to users' experience.

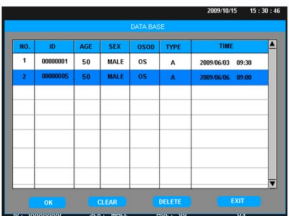


Patient Management

Built-in data archive for permanent memory of 180 patient records.

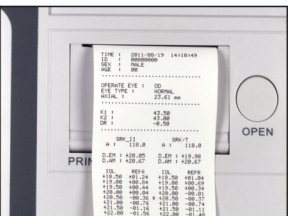



Instant Printout

Instant printout by built-in thermal printer. User-defined printing options.

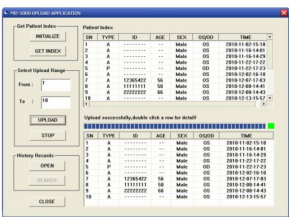


PC Connection

Uploading application software is developed to enable communication with PC and massive storage capability.

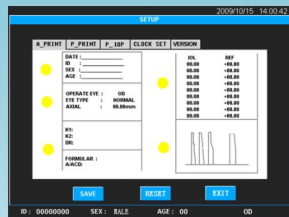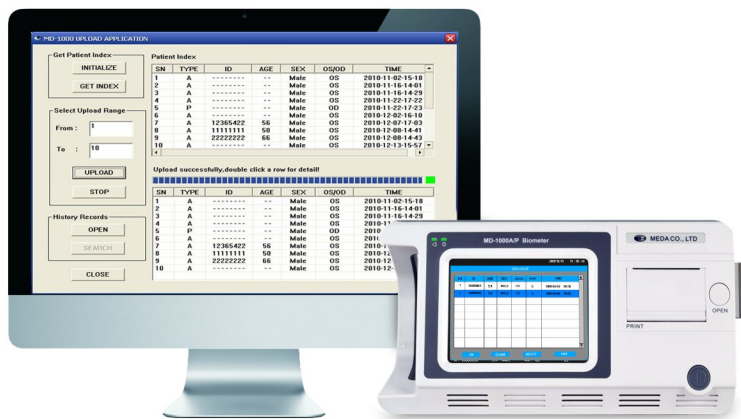


User-defined Interface


User may define acoustic velocity, IOP parameters, printing option etc.

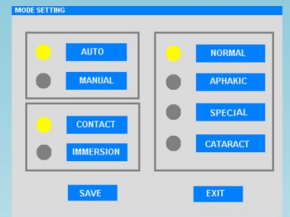


ACCESSORIES

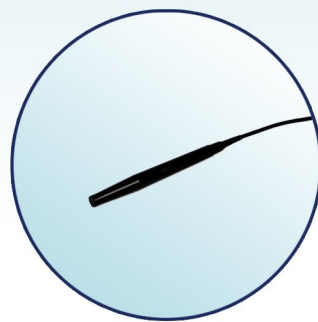

10MHz A-probe

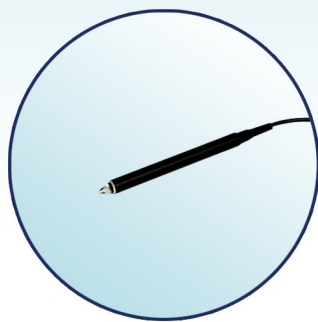


20MHz straight P-probe

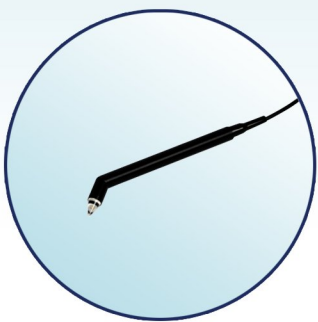

20MHz angled P-probe



Footswitch

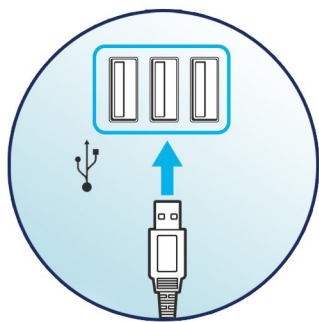

PC Suite